

Data Sheet

pYSG-IBA5

Cat. No.: 5-4605-001

Version: 2.2

Lot No.: 4605-

Revision Date: 04.03.2020

Description	StarGate® Acceptor Vector designed for high-level expression in yeast containing the following elements: <ul style="list-style-type: none"> • Copper inducible promoter (CUP1) for controlled high-level expression • URA3 auxotrophy marker for selection after transformation (do not use URA3 for selection during expression) • LEU2d auxotrophy marker for selection to increase plasmid copy number for expression (do not use LEU2d for selection after transformation) • 2μ ori for episomal replication in yeast • The expressed recombinant protein will be localized in the cytoplasm.
Yeast Expression	Cultivate transformed yeast cells under LEU2d selection until OD600 reaches 0.8 – 1.2 absorbance units. Induce protein expression by addition of copper sulphate to a final concentration of 0.5 mM.
Affinity tag	Strep-Tactin affinity tag (Strep-tag® II) for the purification of recombinant protein via Strep-Tactin resins. The Strep-tag® II is fused to the N-terminus of the recombinant protein.
Resistance	Ampicillin
Form	5 μ g, dissolved in 20 μ l TE buffer, pH 8,0: 10 mM Tris-HCl, 1 mM EDTA
Concentration	250 ng/ μ l
Stability	12 months after shipping
Storage	recommended: 2-8 °C for frequent usage, -20 °C for long-term storage
Shipping	room temperature
Hazards	Product is not classified as hazardous according to (EC) No 1272/2008 [CLP]. A Material Safety Data Sheet is provided.

Note: The sequences have been compiled from information in the sequence database, published literature, and other sources, together with partial sequences obtained by IBA, however, the vectors have not been completely sequenced.



Go digital and help the environment. Please download all up-to-date manuals, protocols and other material from <http://www.iba-lifesciences.com>.

For research use only


Important licensing information

This product is based on StarGate and Strep-tag technologies covered by intellectual property (IP) rights and on completion of the sale IBA grants respective Limited Use Label Licenses to purchaser. IP rights and Limited Use Label Licenses for said technology are further described and identified at <http://www.iba-lifesciences.com/patents.html> or upon inquiry at info@iba-lifesciences.com or at IBA GmbH, Rudolf-Wissell-Str. 28, 37079 Goettingen, Germany. By use of this product the purchaser accepts the terms and conditions of all applicable Limited Use Label Licenses.


Trademark information

The owners of trademarks marked by “®” or “TM” are identified at <http://www.iba-lifesciences.com/patents.html>. Registered names, trademarks, etc. used in this document, even when not specifically marked as such, are not to be considered unprotected by law.

Expression cassette of pYSG-IBA5

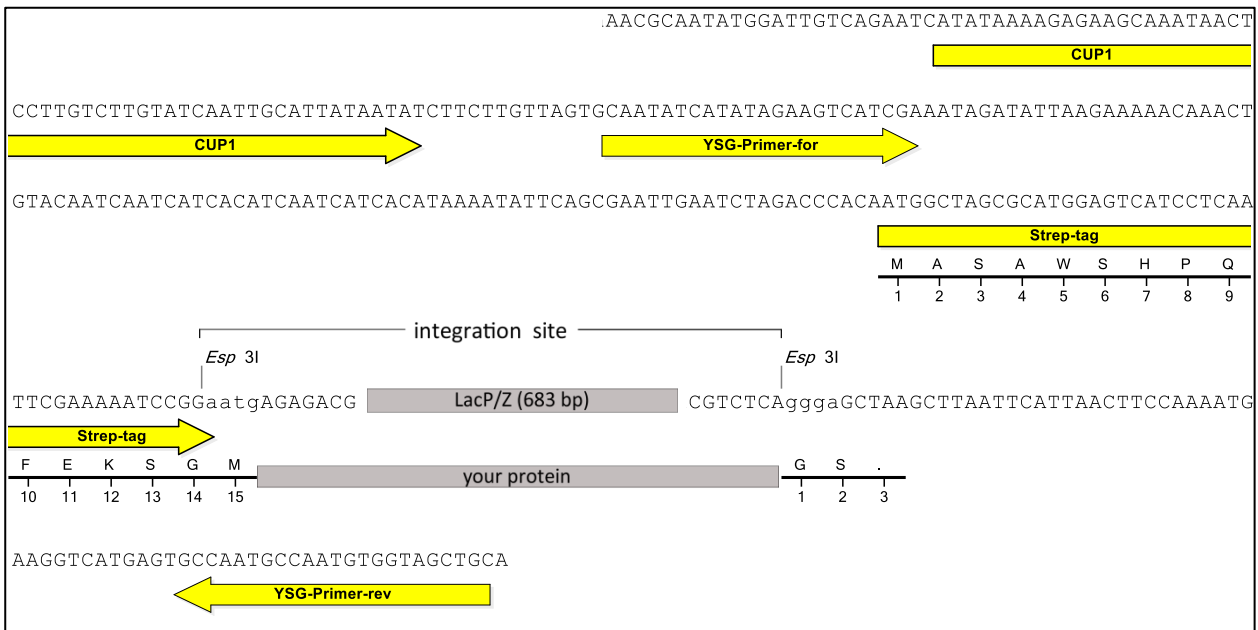

LacP/Z cassette = contains LacZ alpha fragment under control of a separate promoter, which allows alpha complementation of *LacZ* mutations such as *LacZΔM15* as in *E. coli* DH5α or TOP10.
 your protein = after StarGate cloning using *Esp31* your gene of interest will be located here

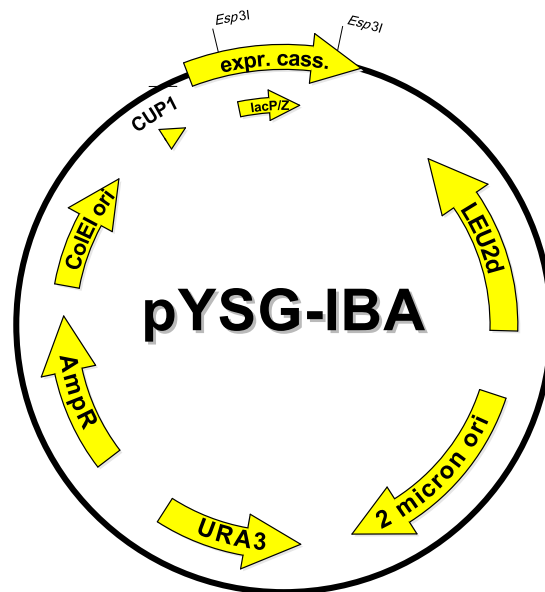

Features	from bp	to bp	Sequencing primer
LEU2d	1668	574	YSG-Primer-for (Cat. No. 5-0000-141) 5' - CAATATCATATAGAAGTCATCGA -3'
2 micron ori	2032	3194	
URA3	4293	3490	YSG-Primer-rev (Cat. No. 5-0000-142) 5' - GCAGCTACCACATTGGCATTGGC -3'
Ampicillin resistance gene	4725	5585	
ColEI ori	5756	6345	
CUP1 promoter	6873	6925	
forward primer binding site	6939	6961	
Strep-tag	7049	7090	
LacZ alpha fragment	7319	7720	
reverse primer binding site	7830	7852	
total vector length		7853	