

IBA Lifesciences GmbH

Rudolf-Wissell-Str. 28 37079 Goettingen Germany

Tel.: +49 (0) 551-5 06 72-0

E-mail: info@iba-lifesciences.com www.iba-lifesciences.com

Data Sheet

pASG-IBA3

Cat. No.: 5-4003-001

Version: 3.0

Revision Date: 14.07.2021

Description	StarGate Acceptor Vector for bacterial expression. The expression cassette is under transcriptional control of the tetracycline promoter/operator. The expressed recombinant protein will be localized in the cytoplasm.		
Affinity tag	Strep-tag® II is fused to the C-terminus of the recombinant protein.		
Cloning Strategy	Cloning into StarGate Acceptor Vectors has to be done with the restriction enzyme Esp3I. There is no Multiple Cloning Site (MCS) available that can be used for the integration of the gene of interest instead (see manual).		
Expression strain	Any E. coli strain. The tet-promoter works independently from the genetic background of E. coli.		
Bacterial Expression	Expression is induced upon addition of 200 μg anhydrotetracycline per 1 liter <i>E. coli</i> shaking culture (A ₅₅₀ = 0.5).		
Resistance	Ampicillin		
Form	5 μg, dissolved in 20 μl TE buffer, pH 8.0: 10 mM Tris/HCl, 1 mM EDTA		
Concentration	250 ng/μl		
Stability	12 months after shipping		
Storage	recommended: 2-8 °C for frequent usage, -20 °C for long-term storage		
Shipping	room temperature		
Hazards	Product is not classified as hazardous according to (EC) No 1272/2008 [CLP]. A Material Safety Data Sheet is provided.		

For research use only


Trademark information

The owners of trademarks marked by """ or "TM" are identified at http://www.iba-lifesciences.com/patents.html. Registered names, trademarks, etc. used in this document, even when not specifically marked as such, are not to be considered unprotected by law.

Important licensing information

This product is covered by intellectual property (IP) rights and on completion of the sale IBA Lifesciences grants respective Limited Use Label Licenses to purchaser. IP rights and Limited Use Label Licenses for said technology are further described and identified at http://www.iba-lifesciences.com/patents.html or upon inquiry at info@iba-lifesciences.com or at IBA Lifesciences GmbH, Rudolf-Wissell-Str. 28, 37079 Goettingen, Germany. By use of this product the purchaser accepts the terms and conditions of all applicable Limited Use Label Licenses.

Expression cassette of pASG-IBA3

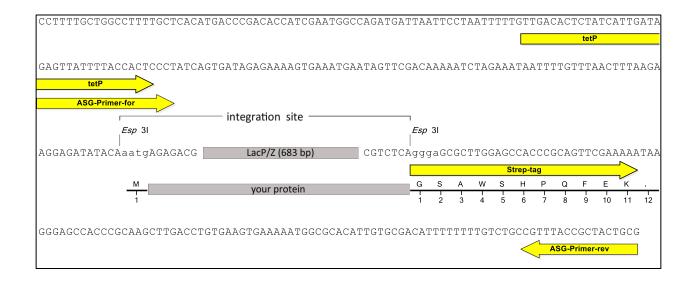


LacP/Z cassette = contains LacZ alpha fragment under control of a separate promoter, which allows alpha

separate promoter, which allows alpha complementation of *LacZ* mutations such as *LacZΔM15* as in *E. coli* DH5α or TOP10.

your protein = after StarGate cloning using Esp3l your gene of

interest will be located here

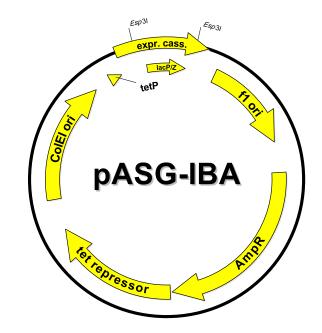

Features	from bp	to bp	Sequencing primer
f1 origin	13	451	ASG-Primer-for
AmpR resistance gene	600	1460	
Tet-repressor	1470	2093	5'- GAGTTATTTTACCACTCCCT -3'
ColEl ori	2246	2834	
Tet promoter	2939	2975	ASG-Primer-rev
forward primer binding site	2959	2978	ASG-Primer-rev
LacZ alpha fragment	3290	3691	5'- CGCAGTAGCGGTAAACG -3'
Strep-tag®II	3755	3787	
reverse primer binding site	3861	3877	
total vector length		3877	